

Kinetik ET-100A UK version Quick Start Guide (Double-sided card)
Size: 50(W) x 75(H) mm [Cover and Back Cover Size: 54(W) x 79(H) mm]
Date: 20201215


Colour : Pantone 173 ■
(Black)

Kinetik Logo Colour : 
Pantone 382C ■ **kinetik**
WELLBEING
(CMYK reference is C30 M0 Y100 K0)

Pantone Cool Grey 11 ■
(CMYK reference is C44 M34 Y22 K77)

50mm

75mm

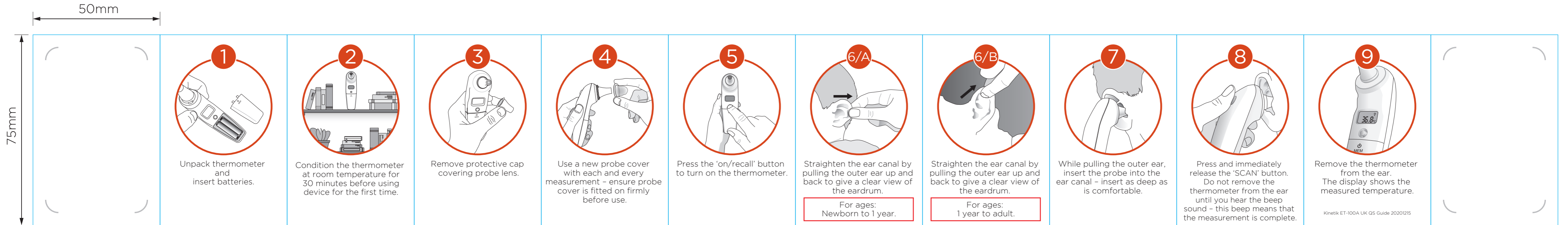

1 Unpack thermometer and insert batteries.

2 Condition the thermometer at room temperature for 30 minutes before using device for the first time.

3 Remove protective cap covering probe lens.

4 Use a new probe cover with each and every measurement - ensure probe cover is fitted on firmly before use.

5 Press the 'on/recall' button to turn on the thermometer.

6/A Straighten the ear canal by pulling the outer ear up and back to give a clear view of the eardrum.
For ages: Newborn to 1 year.

6/B Straighten the ear canal by pulling the outer ear up and back to give a clear view of the eardrum.
For ages: 1 year to adult.

7 While pulling the outer ear, insert the probe into the ear canal - insert as deep as is comfortable.

8 Press and immediately release the 'SCAN' button. Do not remove the thermometer from the ear until you hear the beep sound - this beep means that the measurement is complete.

9 Remove the thermometer from the ear. The display shows the measured temperature.

Kinetik ET-100A UK QS Guide 20201215

54mm

79mm



In association with **St John Ambulance** | **kinetik** WELLBEING

-Quick Start Guide-

ET-100A

Inner Ear Thermometer

Cover

54mm

79mm



kinetik WELLBEING

for more information visit :
kinetikwellbeing.com

Back Cover