

In association with

St John
Ambulance





kinetik
WELLBEING

-Quick Start Guide-

ET-100A

Inner Ear Thermometer

1

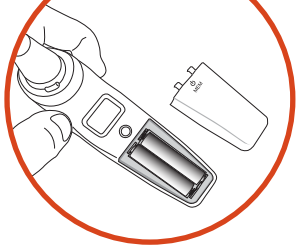

Unpack thermometer
and
insert batteries.

2




Condition the thermometer at room temperature for 30 minutes before using device for the first time.

3




Remove protective cap covering probe lens.

4

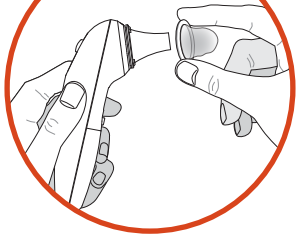


Use a new probe cover with each and every measurement – ensure probe cover is fitted on firmly before use.

5

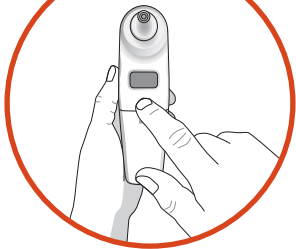


Press the 'on/recall' button to turn on the thermometer.

6/A

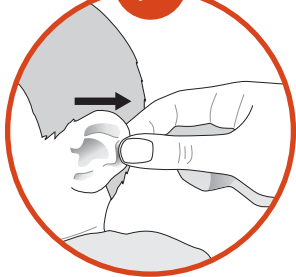

Straighten the ear canal by pulling the outer ear up and back to give a clear view of the eardrum.

For ages:
Newborn to 1 year.

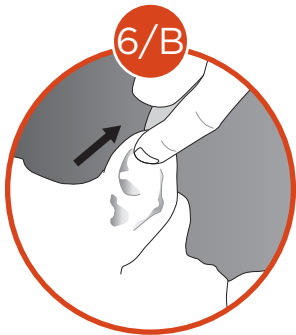

Straighten the ear canal by pulling the outer ear up and back to give a clear view of the eardrum.

For ages:
1 year to adult.



While pulling the outer ear,
insert the probe into the
ear canal – insert as deep as
is comfortable.

8



Press and immediately release the 'SCAN' button.

Do not remove the thermometer from the ear until you hear the beep sound - this beep means that the measurement is complete.



Remove the thermometer
from the ear.
The display shows the
measured temperature.



kinetik
WELLBEING



MIX

Packaging from
sustainable sources

FSC® C166649

for more information visit :
kinetikwellbeing.com